ONE TOUCH

BLOOD PRESSURE MONITOR

UPPER ARM

UA-660CN

Instruction Manual
ORIGINAL

Table of Contents

1. Introduction	E-3
Intended Use / FCC Compliance / Precautions	
About Blood Pressure / About Hypertension	
About Irregular Heartbeat / IHB/AFib Indicator	
Tips for Taking Blood Pressure	
AHA Classification	
2. Getting to Know Your Monitor	E-10
Parts Identification / Display	
Symbols	
3. Initial Setup	E-14
Installing and Changing Batteries	
Connecting Air Hose	
Adjusting Clock	
4. Taking Your Blood Pressure	E-17
Applying and Removing the Cuff	
Normal Measurement	
Measurement with Desired Systolic Pressure	
Recalling and Deleting Memory	
5. Maintenance	E-22
6. Troubleshooting	E-22
7. Optional Accessories	E-24
8. Warranty	E-25
9. Technical Data	E-27

Congratulations on purchasing a state-of-the-art A&D blood pressure monitor. Designed for ease of use and accuracy, this monitor will facilitate your daily blood pressure regimen. We recommend that you read through this manual carefully before using the monitor for the first time.

Intended Use

- The monitor is designed for use on adults only. Do not use on newborns or infants.
- Environment for use: This monitor is intended for home healthcare environment.
- This monitor is designed to measure blood pressure and pulse rate of people for diagnosis.

FCC Compliance Information

This monitor complies with Part 15 of FCC rules. Operation is subject to the following two conditions:

1. This monitor may not cause harmful interference, and
2. This monitor must accept any interference received, including interference that may cause undesired operations.

Precautions

- Precision components are used in the construction of this monitor. Extremes in temperature, humidity, direct sunlight, shock or dust should be avoided.
- Clean the monitor and cuff with a dry, soft cloth or a cloth dampened with water and a neutral detergent. Never use alcohol, benzene, thinner or other harsh chemicals to clean the monitor or cuff.
- Avoid tightly folding the cuff or storing the hose tightly twisted for long periods, as such treatment may shorten the life of the components.
- The monitor and cuff are not water resistant. Prevent rain, sweat and water from soiling the monitor and cuff.

Introduction

- Measurements may be distorted if the monitor is used close to televisions, microwave ovens, cellular telephones, X-ray or other devices with strong electrical fields.
- When using the monitor, confirm that the monitor is clean.
- Used equipment, parts and batteries are not treated as ordinary household waste, and must be disposed of according to the applicable local regulations.
- Do not modify the monitor. It may cause accidents or damage to the monitor.
- To measure blood pressure, the arm must be squeezed by the cuff hard enough to temporarily stop blood flow through the artery. This may cause pain, numbness or a temporary red mark to the arm. This condition will appear especially when measurement is repeated successively. Any pain, numbness, or red marks will disappear with time.
- Take care to avoid accidental strangulation of babies or infants with the hose and cable.
- Do not twist the air hose during measurement. This may cause injury due to continuous cuff pressure.
- Wireless communication devices, such as home networking devices, mobile phones, cordless phones and their base stations, walkie-talkies can affect this blood pressure monitor. Therefore, a minimum distance of 30 cm (12") should be kept from such devices.
- Measuring blood pressure too frequently may cause harm due to blood flow interference. Check that the operation of the monitor does not result in prolonged impairment of blood circulation, when using the monitor repeatedly.
- If you have had a mastectomy or lymph node clearance, please consult a doctor before using the device.
- Do not let children use the monitor by themselves and do not use the monitor in a place within the reach of infants. It may cause accidents or damage.

- There are small parts that may cause a choking hazard if swallowed by mistake by infants.
- Do not apply the cuff on an arm in which another medical device is attached. The equipment may not function properly.
- People who have a severe circulatory deficit in the arm must consult a doctor before using the monitor, to avoid medical problems.
- Do not self-diagnose the measurement results and start treatment by yourself. Always consult your doctor for evaluation of the results and treatment.
- Do not apply the cuff on an arm with an unhealed wound.
- Do not apply the cuff on an arm receiving an intravenous drip or blood transfusion. It may cause injury or accidents.
- Do not use the monitor where flammable gases such as anesthetic gases are present. It may cause an explosion.
- Do not use the monitor in highly concentrated oxygen environments, such as a high-pressure oxygen chamber or an oxygen tent. It may cause a fire or explosion.
- Use of accessories not detailed in this manual may compromise safety.
- Should the battery short-circuit, it may become hot and potentially cause burns.
- Allow the monitor to adapt to the surrounding environment before use (about one hour).
- Clinical testing has not been conducted on newborn infants and pregnant woman. Do not use on newborn infants or pregnant woman.
- Do not touch the batteries and the patient at the same time. That may result in electrical shock.
- Do not inflate without wrapping the cuff around the upper arm.
- When applying the cuff, ensure the arm protector flap is flat against arm to prevent injury to the skin.
- Do not excessively bend or squeeze the air hose.

Introduction

About Blood Pressure

Blood pressure is the force exerted by blood against the walls of the arteries. Systolic pressure occurs when the heart contracts; diastolic pressure occurs when the heart expands. Blood pressure is measured in millimeters of mercury (mmHg).

Blood pressure is affected by many factors: age, weight, time of day, activity level, climate and altitude. Certain activities can alter one's blood pressure. Walking will likely raise blood pressure. Sleeping will likely decrease blood pressure. Not relaxing for several moments before measurements can influence readings as well. In addition to these factors, beverages containing caffeine or alcohol, certain medications, emotional stress and even tight-fitting clothes can make a difference in the readings.

About Hypertension

Hypertension (high blood pressure) is the diagnosis given when readings consistently rise above normal. Hypertension can lead to stroke, heart attack or other illness if left untreated. Referred to as a "silent killer" because it does not always produce symptoms that alert you to the problem, hypertension is treatable when diagnosed early.

In many individuals, hypertension can be controlled by altering lifestyle and minimizing stress, and by appropriate medication prescribed and monitored by your doctor. The American Heart Association recommends the following lifestyle suggestions to prevent or control hypertension:

- Do not smoke
- Reduce salt and fat intake
- Maintain proper weight
- Exercise regularly
- Have regular physical checkups

Blood pressure readings taken in a clinical setting might be elevated as a result of apprehension and anxiety. This is commonly called "white coat hypertension." Self-measurement at home supplements your doctor's readings and provides a more accurate, complete blood pressure history.

In addition, clinical studies have shown that the detection and treatment of hypertension is improved when patients both consult their physicians and monitor their own blood pressure at home.

IHB/AFib Indicator

When the monitor detects an irregular rhythm during the measurements, the IHB/AFib indicator will appear on the display with the measurement values.


Note: We recommend contacting your physician if you see this  IHB/AFib indicator frequently.

What is AFib?

The heart contracts due to electrical signals occurring in heart and sends blood through the body. Atrial fibrillation (AFib) occurs when the electrical signal in the atrium becomes confused and leads to disturbances in the pulse interval. AFib can cause blood to stagnate in the heart, which can easily create clots of blood, a cause of stroke and heart attack.

Pressure Bar Indicator

The indicator monitors the progress of pressure during measurement

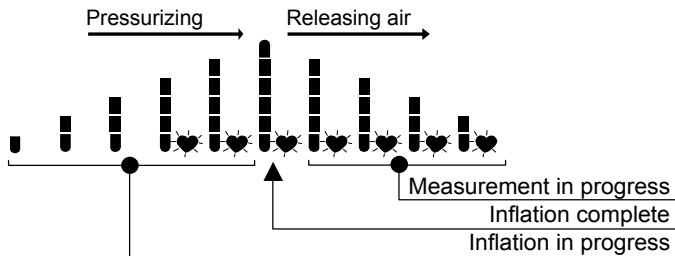Introduction

This blood pressure monitor is designed to detect the pulse and to inflate the cuff to a systolic pressure level automatically.

Tips for Taking Blood Pressure

Before Your Measurement:

- For 30 minutes prior to taking your blood pressure
 - » Do not exercise
 - » Do not drink coffee, caffeinated soda or alcohol
 - » Do not smoke
- Sit quietly for 5-10 minutes
- Ensure you are using the right size cuff

During Your Measurement:

- Do not talk
- Sit with your back straight and supported
- Uncross your legs and place feet flat on floor
- Rest arm on a table so that the cuff is at heart level
- Measure two times a day, in the morning and evening

Notes for Accurate Measurement

- This monitor bases its measurements on the heartbeat. If you have a very weak or irregular heartbeat, the monitor may have difficulty determining your blood pressure.
- Should the monitor detect a condition that is abnormal, it will stop the measurement and display an error symbol. Refer to the section "Symbols" for the description of the symbols.
- If you have emotional stress, the measurement will reflect this stress as a higher (or lower) than normal blood pressure reading and the pulse reading will usually be faster than normal.
- An individual's blood pressure varies constantly, depending on what you are doing and what you have eaten. What you drink can have a very strong and rapid effect on your blood pressure.

AHA Classification

Each segment of the bar indicator corresponds to the AHA blood pressure classification.

Example

Hypertension Stage 2	Hypertension Stage 1	Elevated
		

AHA Classification Indicator


	← Hypertension Stage 2
	← Hypertension Stage 1
	← Elevated
	← Normal

: The indicator displays a segment, based on the current data, corresponding to the AHA classification.


Blood Pressure Classification	Systolic (mmHg)		Diastolic (mmHg)
Normal	< 120	and	< 80
Elevated	120 – 129	and	< 80
Hypertension Stage 1 (High Blood Pressure)	130 – 139	or	80 – 89
Hypertension Stage 2 (High Blood Pressure)	≥ 140	or	≥ 90
HYPERTENSIVE CRISIS (consult your doctor immediately)	≥ 180	and / or	≥ 120

Getting to Know Your Monitor

Parts Identification

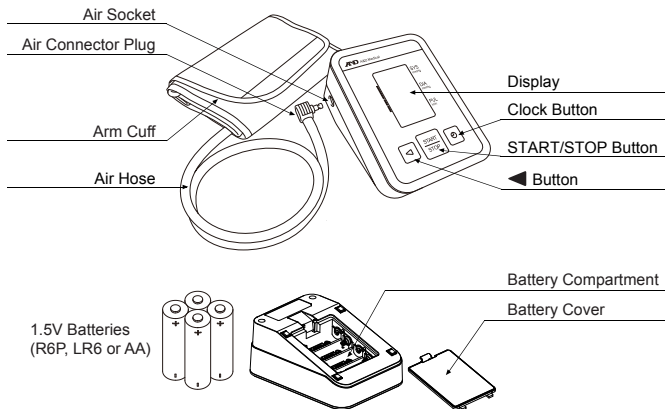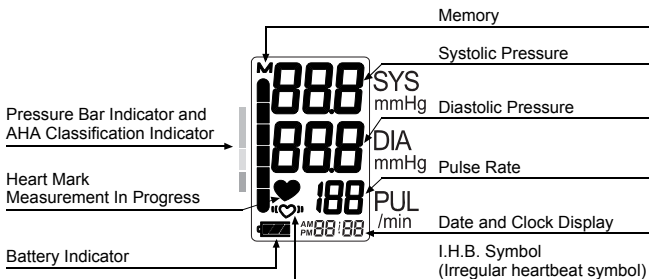Display

Getting to Know Your Monitor

Symbols

Symbols that appear on the display

Symbols	Function/Meaning	Recommended Action
	Appears while measurement is in progress. It blinks when the pulse is detected.	Measurement is in progress. Remain as still as possible.
	Irregular Heartbeat symbol (I.H.B. / AFib) Appears when an irregular heartbeat is detected. It may light when a very slight vibration like shivering or shaking is detected.	_____
M	Previous measurements stored in memory.	_____
	FULL BATTERY – The battery power indicator during measurement.	_____
	LOW BATTERY The battery is low when it blinks.	Replace all batteries with new ones, when the indicator blinks.
AM	Time in the morning	_____
PM	Time in the afternoon	_____

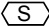
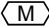
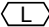
Getting to Know Your Monitor


Symbols printed on the monitor case

Symbols	Function/Meaning
	Standby and turn the monitor on
SYS	Systolic blood pressure in mmHg
DIA	Diastolic blood pressure in mmHg
PUL/min	Pulse per minute
	Battery installation guide
IP	International protection symbol
	Global trade item number
	Lot number
	Serial number
	Refer to instruction manual/booklet
	Keep dry
	Clock setting
	Clock adjustment and memory recall

Getting to Know Your Monitor

Symbols that are printed on the cuff


Symbols	Function/Meaning	Recommended Action
●	Artery position mark	Set the ● mark on the artery of the upper arm or in line with the ring finger on the inside of the arm.
▲	Index	_____
REF	Catalog number	_____
LOT	Lot number	_____
	Type BF: applied part	_____
W	Proper fit range for the Wide Range cuff.	_____
L	Proper fit range for the Large cuff.	_____
M	Proper fit range for the Medium cuff.	_____
S	Proper fit range for the Small cuff.	_____
	Range to use the Small cuff. Under range printed on the Medium cuff and Wide Range cuff.	Use the Small cuff
	Range to use the Medium cuff. Under range printed on the Large cuff and over range printed on the Small cuff.	Use the Medium cuff
	Range to use the Large cuff. Over range printed on the Medium cuff and Wide Range cuff.	Use the Large cuff

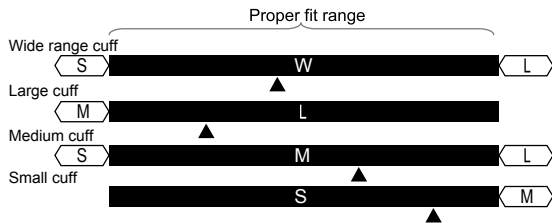Initial Setup

Installing/Changing the batteries

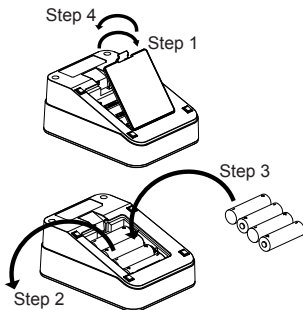1. Remove the battery cover.
2. Remove the used batteries from the battery compartment when changing them.
3. Insert new batteries or replace used batteries as shown, taking care that the polarities (+) and (-) are correct.
4. Replace the battery cover. Use only R6P, LR6 or AA batteries.


CAUTION

- Insert the batteries as shown in the battery compartment. If installed incorrectly, the monitor will not work.
- When  (LOW BATTERY mark) blinks on the display and the monitor announces that the battery needs to be replaced, replace all batteries with new ones. Do not mix old and new batteries. It may shorten the battery life, or cause the monitor to malfunction.
- Replace the batteries two seconds or more after the monitor turns off.
- If  (LOW BATTERY mark) does not appear then the batteries are drained.
- The battery life varies with the ambient temperature and may be shorter at low temperatures.
- Remove the batteries if the monitor is not to be used for a long time. The batteries may leak and cause a malfunction.
- Use the specified batteries only. The batteries provided with the monitor are for testing the monitor performance and may have limited life.

Connecting the Air Hose

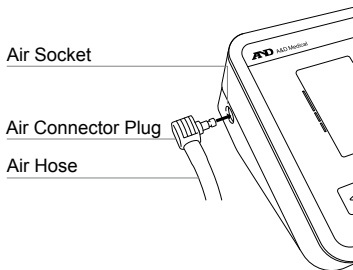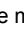
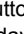
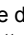
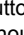
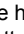
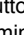
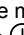
Insert the air connector plug into the air socket firmly.

Initial Setup


Adjusting the build-in clock

Adjust the clock prior to use.

1. Press the , the year starts blinking.
2. Select the year using the  button. Press the  button to set the current year and move to month/day selection. The date can be set anywhere between the years 2024 and 2079.
3. Select the month using the  button. Press the  button to set the current month and move to day selection.
4. Select the day using the  button. Press the  button to set the current day and move to hour/minute selection.
5. Select the hour using the  button. Press the  button to set the current hour and move to minute selection.
6. Select the minute using the  button. Press the  button to turn the monitor off.

Note: After three minutes of non-operation, the monitor will turn off automatically. When the clock has not been set, the clock display indicates dashes as shown to the right.

When using the monitor for the first time, the clock is not adjusted. When the monitor is disconnected from the power supply, the set date and time will be erased. When the set date and time is erased, please adjust again.

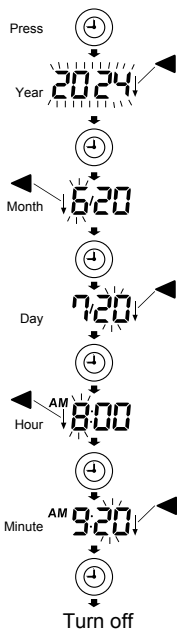Pressing the **START/STOP** button will turn the monitor off anytime

-/--

-:--



Taking Your Blood Pressure

Applying the arm cuff

1. Insert arm (preferably the left arm) into cuff as shown in the figure to the right.
2. Wrap the cuff around the upper arm, about 0.5" (1-2 cm) above the inside of the elbow, as shown in the figure to the right.
3. Place the cuff directly against the skin, as clothing may cause a faint pulse and result in a measurement error. Constriction of the upper arm, caused by rolling up a shirtsleeve, may prevent accurate readings.
4. Confirm the arm protector is flat against the arm as shown in the figure below.
5. Confirm that the index ▲ points within the proper fit range.


Do not roll up shirtsleeve tightly

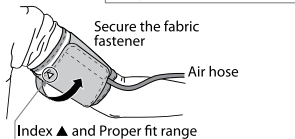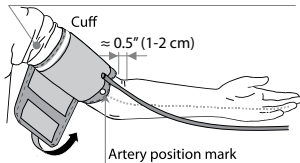
Note: During measurement, it is normal for the cuff to feel very tight. (Do not be alarmed).

Taking Your Blood Pressure

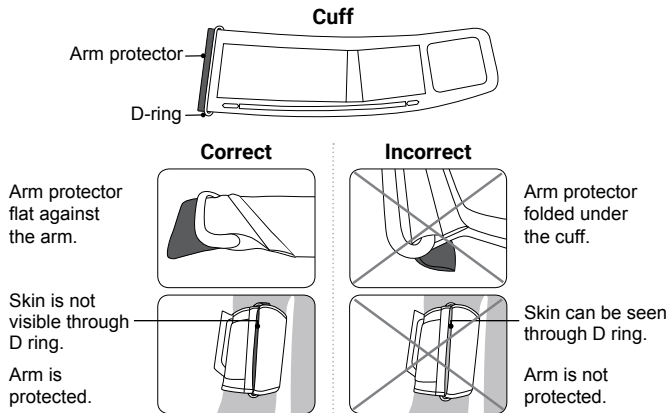Removing the arm cuff

1. Wait for the cuff to deflate
2. Unfasten the Velcro

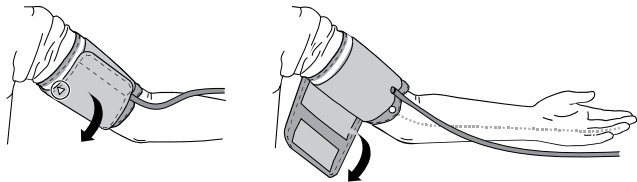3. Remove the cuff

Taking Your Blood Pressure

Normal Measurement

1. Place the cuff on the arm (preferably the left arm). Sit quietly during measurement.
2. Press the START/STOP button. All of the display segments are displayed. Zero is displayed blinking briefly. Then the display changes, as indicated in the figure at the right, as the measurement begins. The cuff starts to inflate. It is normal for the cuff to feel very tight. A pressure bar indicator is displayed, as in the figure at the right, during inflation.


Note: If you wish to stop inflation at any time, press the START/STOP button again.

3. When inflation is complete, deflation starts automatically and the ♥ (heart mark) blinks, indicating that the measurement is in progress. Once the pulse is detected, the mark blinks with each pulse beat.

Note: If an appropriate pressure is not obtained, the monitor starts to inflate again automatically.

4. When the measurement is complete, the systolic and diastolic pressure readings and pulse rate are displayed. The cuff exhausts the remaining air and deflates completely.
5. Press the START/STOP button again to turn off the power.

Note: This product is provided with an automatic power shut-off function, which turns the power off approximately one minute after measurement.

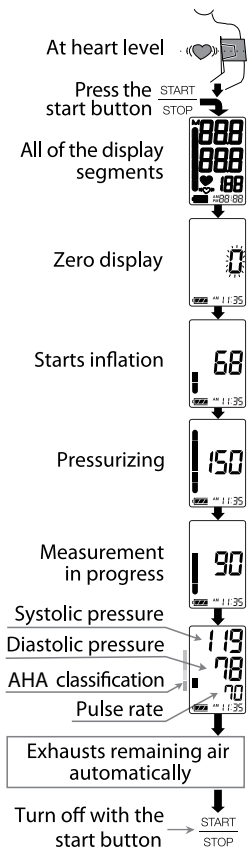
Taking Your Blood Pressure

Measurement with the Desired Systolic Pressure

If your systolic pressure is expected to exceed 230 mmHg, use this procedure.

Use this method when re-inflation occurs repeatedly or when the results are not displayed even if the pressure decreases to 20 mmHg or less.

1. Place the cuff on the arm (preferably the left arm).
2. Press and hold the START/STOP button until a number about 30 to 40 mmHg higher than your expected systolic pressure appears.
3. When the desired number is reached, release the START/STOP button to start measurement. Continue to measure your blood pressure as described on the previous page.

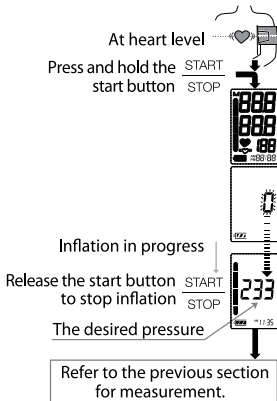
Taking Your Blood Pressure

Recalling the memory

Note: This monitor stores the last 60 measurements in memory.

1. Press the ◀ button. The average of the last three measurements is displayed. (If no data, "0" is displayed. Press the ◀ or START/STOP button to turn the device off.)



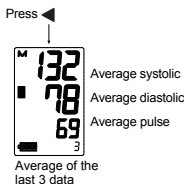
2. Each time the ◀ button is pressed, the memory data is displayed as follows. Most recent data (No.n, in the example, No.35)

Three seconds after the data number display, the measurement data is displayed.


Last data (No.1). Three seconds after the data number display, the measurement data is displayed.

3. After the last data is displayed, press the ◀ button to return to the average of the last 3 measurements.
4. Press the START/STOP button to turn the device off. After one minute of non-operation, the device will turn off automatically.

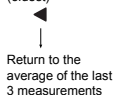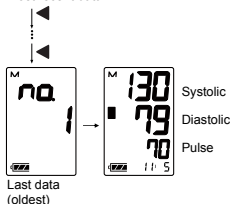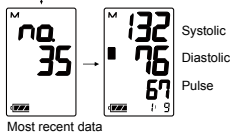Deleting Data Stored in Memory

When the monitor is turned off, press and hold the ◀ button. The monitor displays the M mark, deletes data stored in memory when the M mark is blinking and turns off automatically.

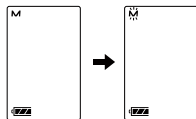Maintenance

Do not open the monitor. It uses delicate electrical components and an intricate air unit that could be damaged. If you cannot fix the problem using the troubleshooting instructions, contact the authorized dealer in your area or our customer service department. A&D customer service can provide technical assistance and spare parts.

Troubleshooting

Problem	Possible Reason	Recommended Action
Nothing appears in the display, even when the power is turned on.	Batteries are drained.	Replace all batteries with new ones.
	Battery terminals are not in the correct position.	Reinstall the batteries with negative and positive terminals matching those indicated on the battery compartment.
The cuff does not inflate.	Battery voltage is too low.  (LOW BATTERY mark) blinks. If the batteries are drained completely, the mark does not appear.	Replace all batteries with new ones.
The monitor does not measure. Readings are too high or too low.	The cuff is not applied properly.	Apply the cuff correctly.
	You moved your arm or body during measurement.	Make sure you remain very still and quiet during measurement.
	The cuff position is not correct.	Sit comfortably and still. Place your arm on a table with your palm facing upward and the cuff at the same level as your heart.
	_____	If you have a very weak or irregular heartbeat, the monitor may have difficulty in determining your blood pressure.

Troubleshooting

Problem	Possible Reason	Recommended Action
Other	The value is different from that measured at a clinic or doctor's office.	At a clinic or doctor's office apprehension may cause an elevated reading. Home measurement reduces the effects of outside influences on blood pressure readings, supplementing the doctor's readings.
	_____	Remove the batteries. Place them back properly and try the measurement again.

Note: If the actions described above do not solve the problem, contact the dealer. Do not attempt to open or repair this product, as any attempt to do so will make your warranty invalid.

Symbols	Function/Meaning	Recommended Action
Err	Unstable blood pressure due to movement during the measurement.	Take another measurement. Remain still during the measurement.
	The systolic and diastolic values are within 10 mmHg of each other.	
	The pressure value did not increase during inflation.	
Err UF	The cuff is not applied correctly.	Apply the cuff correctly, and try the measurement again.
E	Pulse display error. The pulse is not detected correctly.	
Err E	Monitor internal error	Press the START/STOP button. If the error still appears, contact the dealer.
Err 9		

Optional Accessories

Accessories are sold separately

Part Number	Description	Specification
UA-420A	Wide Range Cuff	8.6-16.5" (22-42 cm)
UA-291A	Large Cuff	12.2-17.7" (31-45 cm)
UA-290A	Medium Cuff	9.0-14.6" (23-37 cm)
UA-289A	Small Cuff	6.3-9.4" (16-24 cm)
UA-372A	Preformed Wide Range Cuff	8.6-16.5" (22-42 cm)
<i>Arm size: The circumference at the biceps.</i>		

LIMITED WARRANTY

A&D Medical

For purchasers within the US and Canada only:

Product	Consumer Warranty Term
Monitor UA-660CN	5 Years
Cuff	2 Years

For purchases made outside the United States or Canada, please contact the local distributor or retailer.

Limited Warranty:

A&D Medical (“A&D”) warrants to the first purchaser (“You”) that the A&D product You purchased (the “Product”) will be free from defects in material, workmanship and design for the applicable Warranty Term stated above from the date You purchased the Product under normal use. This Limited Warranty is personal to You and is not transferable. If the Product is defective, then You return the Product to A&D in accordance with the procedure set forth below. A&D’s warranty obligation is limited to the repair or replacement, at A&D’s option, of the defective Product that has been returned by You within the warranty period. Such repair or replacement will be at no charge to You. The repaired or replacement Product is warranted here-under for the longer of the remainder of the original warranty period or 90 days from the date of shipment of the repaired or replacement Product.

To obtain a warranty service, please contact us in **US at 1-888-726-9966** or in **Canada at 1-800-461-0991** for return address, shipping and handling fee, and other instructions for processing warranty. Please ensure you have satisfactory proof of the date of Your purchase and a description of the defect. Returns will not be accepted unless a Return Material Authorization (RMA) Number has been issued from A&D Customer Service Representative.

This Limited Warranty does not cover, and A&D will not be liable for (i) any shipment damage, (ii) any damage or defect due to misuse, abuse, failure to use reasonable care, failure to follow written instructions enclosed with the Product, accident, subjecting the Product

Warranty

conditions, or modification, alteration or repair by anyone other than A&D or persons authorized by A&D, or (iii) expendable or consumable components.

THIS LIMITED WARRANTY IS THE ONLY WARRANTY PROVIDED BY A&D; THERE ARE NO OTHER EXPRESS WARRANTIES. If A&D cannot reasonably repair or replace the Product, A&D will refund the amount You paid for the Product (not including taxes), less a reasonable charge for usage. To receive a refund you must have returned the Product and all associated materials to A&D.

The above remedy of repair, replacement or refund is your only and exclusive remedy. IN NO EVENT SHALL A&D BE LIABLE FOR ANY DAMAGES, INCLUDING, WITHOUT LIMITATION, LOST PROFITS, LOST INFORMATION OR REPLACEMENT COSTS, ARISING OUT OF YOUR USE OF OR INABILITY TO USE THE PRODUCT, INCLUDING, WITHOUT LIMITATION, ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES, EVEN IF A&D HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. Some states do not allow the exclusion of incidental or consequential damages, so that the above exclusions may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that may vary from state to state.

No distributor, dealer or other party is authorized to make any warranty on behalf of A&D or to modify this warranty, or to assume for A&D any liability with respect to its products.

Technical Data

Type	UA-660CN
Measurement method	Oscillometric measurement
Memory	Last 60 measurements
Measurement range	Pressure: 0-299 mmHg Systolic pressure: 60-279 mmHg Diastolic pressure: 40-200 mmHg Pulse: 40-180 beats /min
Static accuracy	Pressure: ± 3 mmHg Pulse: ± 5 %
Rating	DC6V 3W
Power supply	4 x 1.5V batteries (R6P, LR6 or AA)
Number of measurements	Approx. 700 times LR6 (<i>alkaline batteries</i>) Approx. 200 times R6P (<i>manganese batteries</i>) With pressure value 180 mmHg, room temperature 23 °C
Classification	Internally powered ME equipment (<i>Supplied by batteries</i>). Continuous operation mode
Clinical test	ISO81060-2:2020
EMD	IEC 60601-1-2: 2014+A1:2020
Operating conditions	50 to 104 °F / 15 % to 85 %RH +10 to +40 °C / 15 % to 85 %RH 800 to 1060hPa
Transport/Storage conditions	-4 to 140 °F / 10 % to 95 %RH -20 to +60 °C / 10 % to 95 %RH 700 to 1060hPa
Dimensions	Approx. 3.8" (W) x 2.6" (H) x 5.4" (D) Approx. 96 (W) x 67 (H) x 138 (D) mm

Technical Data

Weight	Approx. 9.2 oz., excluding the batteries Approx. 260 g, excluding the batteries
Ingress protection	Device: IP20
Applied part	Cuff Type BF 
Useful life	Monitor: 5 years (<i>when used six times a day</i>) Cuff: 2 years (<i>when used six times a day</i>)

Notes: Specifications are subject to change without prior notice. IP classification is the degrees of protection provided by enclosures in accordance with IEC 60529. This device is protected against solid foreign objects of 12 mm diameter and greater such as a finger. This device is not protected against water.

EMD Technical Data

Battery-operated Blood Pressure Monitor

Medical Electrical Equipment needs special precautions regarding EMD and needs to be installed and put into service according to the EMD information provided in the following.

Portable and mobile RF communication equipment (e.g. cell phones) can affect Medical Electrical Equipment.

The use of accessories and cables other than those specified may result in increased emissions or decreased immunity of the unit.

Table 1 - EMISSION Limits -

Phenomenon	Compliance
Conducted and radiated RF EMISSION CISPR 11	Group 1, Class B

Table 2 - IMMUNITY TEST LEVELS : Enclosure Port -

Phenomenon	IMMUNITY TEST LEVELS
Electrostatic discharge IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Radiated RF EM fields IEC 61000-4-3	10 V/m 80 MHz - 2.7 GHz 80 % AM at 1 kHz
Proximity fields from RF wireless communications equipment IEC 61000-4-3	See table 3
Rated power frequency magnetic fields IEC 61000-4-8	30 A/m 50 Hz or 60 Hz
proximity magnetic fields IEC 61000-4-39	See table 4

Technical Data

Table 3 - Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment -

Test frequency (MHz)	Band (MHz)	Service	Modulation	IMMUNITY TEST LEVEL (V/m)
385	380 - 390	TETRA 400	Pulse modulation 18 Hz	27
450	430 - 470	GMRS 460 FRS 460	FM ± 5 kHz deviation 1 kHz sine	28
710	704 - 787	LTE Band 13, 17	Pulse modulation 217 Hz	9
745				
780				
810	800 - 960	GSM 800/900 TETRA 800 iDEN 820 CDMA 850 LTE Band 5	Pulse modulation 18 Hz	28
870				
930				
1720	1700 - 1990	GSM 1800 CDMA 1900 GSM 1900 DECT LTE Band 1, 3, 4, 25 UMTS	Pulse modulation 217 Hz	28
1845				
1970				
2450	2400 - 2570	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE Band 7	Pulse modulation 217 Hz	28
5240	5100 - 5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	9
5500				
5785				

Technical Data

Table 4 - Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields-

Test frequency	Modulation	IMMUNITY TEST LEVEL(A/m)
30 kHz	CW	8
134.2 kHz	Pulse modulation 2.1 kHz	65
13.56 MHz	Pulse modulation 50 kHz	7.5

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© 2025 All trademarks are property of their rightful owners
A&D Company, Limited
A&D Engineering, Inc.