

EasyOne® Air

Simple. Certain. Proven.



Innovators for Optimal Lung Care

EasyOne Air is the new, all-in-one revolution in lung function testing. Built for healthcare providers large and small, EasyOne Air puts the power of the most advanced spirometry testing tools in the palm of your hand.

Revolutionary TrueFlow ultrasound technology means unprecedented accuracy for reliable diagnosis and trusted treatment. Equipped with ready-to-use EMR connectivity and a user-friendly interface, the EasyOne Air makes the best in lung care available anywhere.

Learn more about EasyOne Air at nddmed.com

300 Brickstone Square, Suite 604 Andover, MA 01810 | (978) 470-0923 | customerservice@nddmed.com



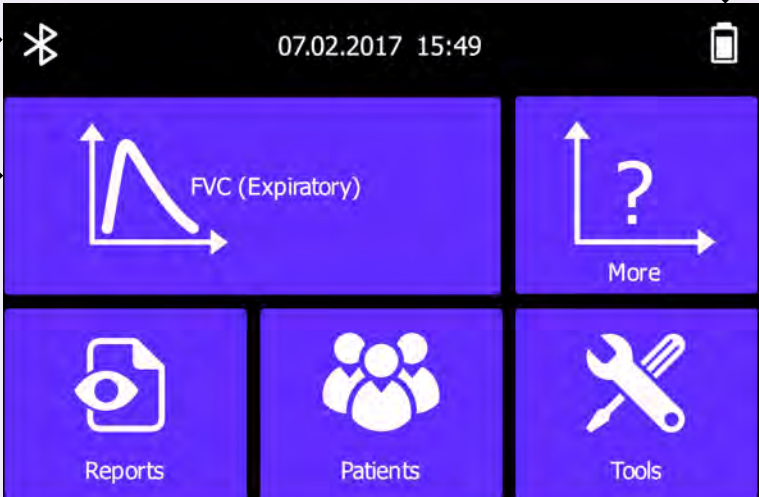
Simple. Certain. Proven.

Rechargeable battery

Wireless connectivity to the pc via Bluetooth



Easy navigation

- Quickly assess test quality with automated QC grades and messaging



- Instant interpretation
- EasyOne Connect for seamless EMR integration
 - Internal storage for up to 10,000 tests
- 3 year exchange warranty
- Calibration free
- Immediate test quality feedback

Large color touch Screen



Incentive screens available with EasyOne Connect

TrueFlow
makes the difference

The original ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity. ndd TrueFlow is a hygienic and resistance-free solution that does not require calibration during its lifetime.

EasyOne Connect
intelligent workflows

EasyOne Connect software offers one database and one platform for all EasyOne point-of-care-solutions. Data management has never been easier.

Specifications

Test Types FVC, FVL, SVC, MVV and Pre/Post Testing

Languages English, French, German, Spanish

Printing Options Direct to printer or with EasyOne Connect software

EMR connectivity EasyOne Connect

Dimensions 3.4X6.1X1.4", 13 OZ

Meets ATS/ERS 2005, OSHA and SSD requirements

Part Numbers

2500-2A EasyOne Air spirometry system

5050-50 EasyOne Flow Tube box of 50

5050-200 EasyOne Flow Tube box of 200