

**A premium Bluetooth® blood pressure monitor with advanced features and voice guidance, making it ideal for a wide range of users, including those who are visually impaired.**

**Product Benefits**

- Audible instructions and results in English, Spanish and French for enhanced accessibility
- TriCheck™ minimizes the impact of individual measurement fluctuations by automating 3 sequential measurements and averaging in accordance with clinical guidelines
- Experience accurate readings with advanced cuff fit detection, avoiding errors from incorrect cuff placement
- Ensure precise measurements with integrated movement detection, minimizing the risk of inaccurate results due to movement
- Automatically records and tracks measurements with the highly rated Heart Track app to enable trending and sharing for better health management

**Comfort and Convenience**

- Advanced inflation and deflation technology provides faster and more comfortable readings
- Includes a carrying case and batteries for immediate use and portability

**The A&D Promise**

- Clinical Grade Accuracy delivers measurements that are as precise as those taken in a medical setting
- Extended 5-year warranty on monitor provides long-term peace of mind and greater value



Advanced Averaging

90 Measurement Memory


Irregular Heartbeat Detection (IHB)

Movement Error

Pressure Rating Indicator

Battery Life Indicator

Cuff Fit Error

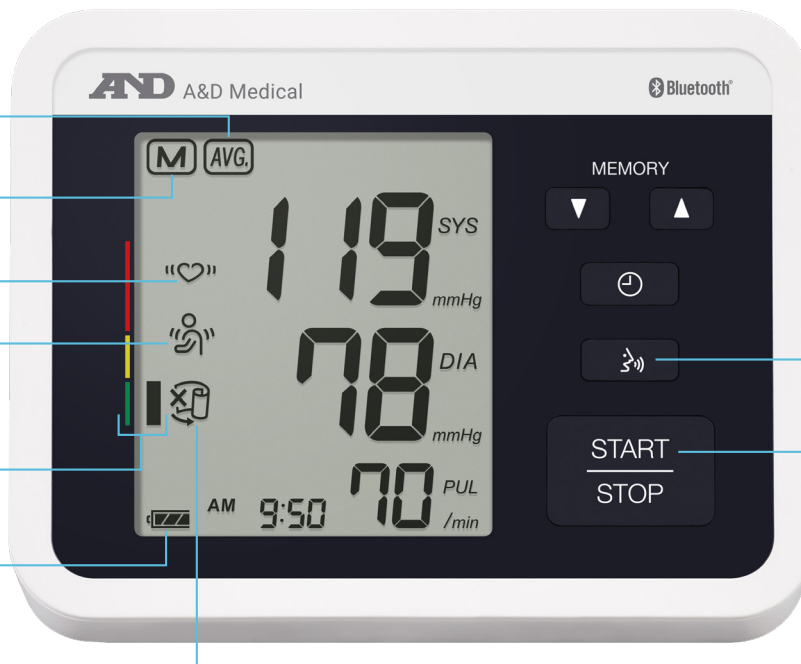

Voice Guided

One Button Operation



**Unlock the Potential**  
Download the FREE  
**A&D Heart Track**  
Companion App



## Technical Specifications

Measurement Method	Oscillometric
Memory	90 measurements
Cuff Size	8.6-16.5" (22-42 cm)
Measurement Range	<b>Pressure:</b> 0-299 mmHg, <b>Systolic pressure:</b> 60-279 mmHg <b>Diastolic pressure:</b> 40-200 mmHg, <b>Pulse:</b> 40-180 beats/min
Measurement Accuracy	<b>Pressure:</b> ±3 mmHg, <b>Pulse:</b> ±5%
Power Supply	4 x 1.5V batteries (R6P, LR6 or AA) or AC adapter (TM-233)
Battery Life	Approx. 700 times LR6 (alkaline batteries) Approx. 200 times R6P (manganese batteries)
Monitor Dimensions	5.5 x 4.6 x 2.5 in (14 x 11.8 x 6.4 cm)
Monitor Weight	11.3 oz (320g) without batteries
Display Size (WxH)	2.4 x 2.8" (6.2 x 7.2cm)
Clinical Test	ISO81060-2: 2020
Operating Conditions	+50°F to +104°F (10°C to 40°C) / 15% to 85% RH / 800 to 1060 hPa
Storage/Transport Conditions	-4°F to +140°F (-20°C to 60°C) / 10% to 95% RH / 700 to 1060 hPa

## Packaging and Shipping Specifications

	Dimensions (WxLxH)	Weight	Case Cube	UPC / GTIN
Individual Unit Box	5.1 x 4.1 x 7.7 in (13 x 10.5 x 19.5 cm)	1.7 lb (770 g)	0.09 ft <sup>3</sup> (0.003 m <sup>3</sup> )	0-93764-60451-5
Case (10 units)	11 x 21.7 x 8.7 in (28 x 55 x 22 cm)	18.5 lb (8.4 kg)	1.2 ft <sup>3</sup> (0.03 m <sup>3</sup> )	1-00-93764-60451-2

## Optional Accessories

Name	Order Code
AC Power Adapter	TB-233

### Compatible Replacement Cuffs- Packaged in individual box

Small Cuff 6.3-9.4" (16-24 cm)	UA-289A
Medium Cuff 9.0-14.6" (23-37 cm)	UA-290A
Large Cuff 12.2-17.7" (31-45 cm)	UA-291A
Wide Range Cuff 8.6-16.5" (22-42 cm)	UA-420A
Cuffs are not made with natural rubber latex	

## What is In The Box

- 1 Blood pressure monitor
- 1 Blood pressure cuff 8.6-16.5" (22-42 cm)
- 1 Instruction manual
- 1 Quick start guide
- 1 Carrying case
- 4 AA batteries