

A wireless blood pressure monitor that automatically records and tracks measurements via mobile app to enable trending and sharing.

Product Benefits

- Displays the average of all measurements providing a quick snapshot in Heart Track mobile app
- Automatically records and tracks 30 measurements with the highly rated Heart Track app to enable trending and sharing for better health management
- Irregular Heartbeat (IHB) indicates the presence of an irregular heartbeat

Comfort and Convenience

- One button operation for quick and easy functionality
- Advanced inflation and deflation technology provides faster and more comfortable readings
- Includes a carrying case and batteries for immediate use and portability

The A&D Promise


- Clinical Grade Accuracy delivers measurements that are as precise as those taken in a medical setting
- Extended 5-year warranty on monitor provides long-term peace of mind and greater value



Bluetooth®

Pressure Rating Indicator

Battery Life Indicator

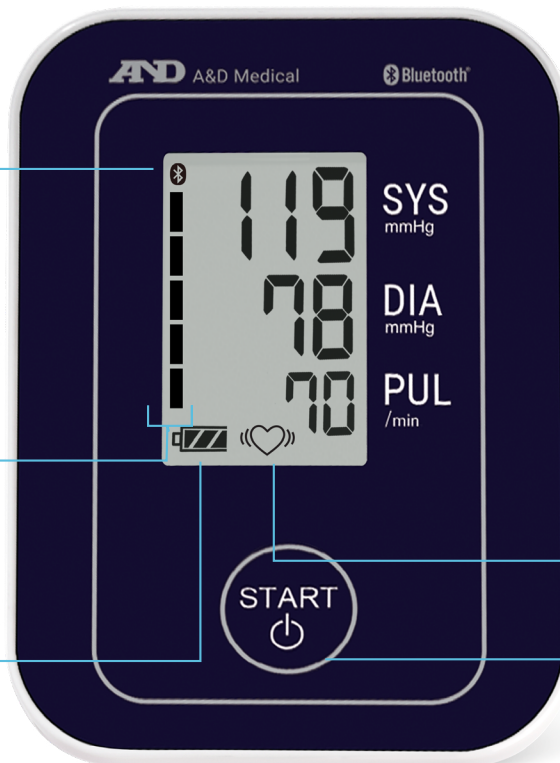

Hypertension
CANADA



Unlock the Potential
Download the FREE
A&D Heart Track
Companion App



Irregular heartbeat
symbol (IHB)/AFIB

One Button
Operation

Technical Specifications

Measurement Method	Oscillometric
Memory	Unlimited on App
Cuff Size	22-42 cm (8.6-16.5")
Measurement Range	Pressure: 0-299 mmHg, Systolic pressure: 60-279 mmHg Diastolic pressure: 40-200 mmHg, Pulse: 40-180 beats / minute
Measurement Accuracy	Pressure: ± 3 mmHg, Pulse: $\pm 5\%$
Power Supply	4 x 1.5V batteries (R6P, LR6 or AA) or AC adapter (TB-233)
Battery Life	Approx. 700 times LR6 (alkaline batteries) Approx. 200 times R6P (manganese batteries)
Monitor Dimensions	9.6 x 13 x 6.8 cm (3.8 x 5.1 x 2.7")
Monitor Weight	249 g (8.8 oz) without batteries
Display Size	3.6 x 5.3 cm (1.4 x 2.1")
Clinical Test	ISO81060-2: 2013
Operating Conditions	+10°C to +40°C (50°F to 104°F) / 15% to 85% RH / 800 to 1060 hPa
Storage/Transport Conditions	-20°C to +60°C (-4°F to 140°F) / 10% to 95% RH / 700 to 1060 hPa

Packaging and Shipping Specifications

	Dimensions (WxLxH)	Weight	Individual Item Cube	UPC / GTIN
UA-651CNBLE Individual Unit (Box)	13 x 10 x 19 cm (5.1 x 4 x 7.7")	0.7 kg (1.52 lb)	0.0024 m ³ (0.09 ft ³)	0-93761-60348-8
UA-651CNBLE Case (10 units)	55 x 28 x 22 cm (21.7 x 11 x 8.5")	7.2 kg (15.87 lb)	0.033 m ³ (1.18 ft ³)	1-00-93764-60348-5

Optional Accessories

Name	Order Code
AC Power Adapter	TB-233

Compatible Replacement Cuffs- Packaged in individual box

Small Cuff 16-24 cm (6.3-9.4")	UA-289A
Medium Cuff 23-37 cm (9.0-14.6")	UA-290A
Large Cuff 31-45 cm (12.2-17.7")	UA-291A
Wide Range Cuff 22-42 cm (8.6-16.5")	UA-420A

Wholesaler

	Codes
McKesson	082959
Kohl & Frisch	159743
LPG	35250
Imperial	490725

Included In The Box

- 1 Blood pressure monitor
- 1 Blood pressure cuff (22-42 cm)
- 1 Instruction manual
- 1 Quick start guide
- 1 Carrying case
- 4 AA batteries

Health Canada License - 94858



6695 Millcreek Drive, Unit 6
Mississauga, ON L5N 5R8, Canada
1-800-461-0991
andmedical.com