

A&D Medical's tabletop blood pressure kiosk is a compact powerhouse ideal for pharmacies, mobile clinics, and workplace screenings. Designed for durability and frequent use, it allows a wide range of patients to easily self-measure their blood pressure and obtain accurate results.


## Where Simplicity Meets Precision

### Flexibility

- Customizable print out for the patient
- Measurements can be uploaded to the highly rated Heart Track app via printable QR code
- Small footprint is ideal for integration in a health screening station or tight spaces
- Extended cuff range fitting arms from 7.1-16.5 in (18-42 cm)
- Optional connectivity configurations

### Simplicity

- One button operation for quick and simple measurements
- No installation required, just remove from box, plug in and turn on power
- Provides a sanitary solution with an easily replaceable anti-bacterial cuff cloth and optional disposable cuff covers

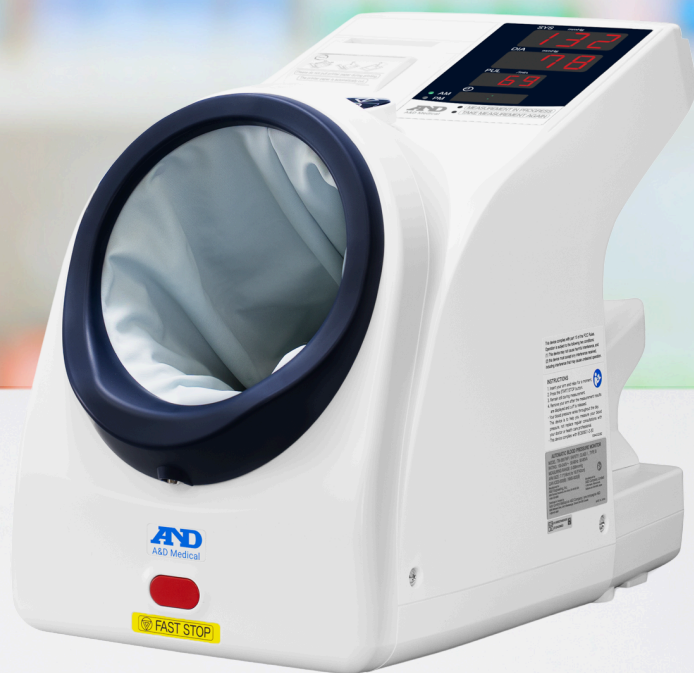

### Unmatched Reliability

With industry-leading reliability and no regular maintenance required, we ensure consistent, trustworthy measurements.



### Clinical Grade Accuracy

Delivers precise measurements that have been validated in a clinical setting and published in a peer reviewed journal.



### Trusted Experts

A&D Medical has nearly 50 years of leadership in developing precision health monitoring technology.

## Technical Specifications

|                              |   |
|------------------------------|---|
| Measurement Method           | Oscillometric   |
| Memory                       | 200 memory  |
| Cuff Size                    | 7.1-16.5 in / 18-42 cm  |
| Display Method               | 3-digit display LED & LED lamp  |
| Time / Date                  | 12 or 24 hour format  |
| Connectivity                 | Optional RS-232C (TM-2657WP-RS)<br>Optional Bluetooth® (TM-2657WP-BT)<br>Optional Bluetooth® Low Energy (TM-2657WP-BLE) |
| Measurement Range            | Pressure: 0-299 mmHg<br>Systolic pressure: 40-279 mmHg<br>Diastolic pressure: 20-200 mmHg<br>Pulse: 30-240 beats/min    |
| Measurement Accuracy         | Static Pressure: ±3 mmHg<br>Pulse: ±5%  |
| AC Power Supply              | 100-240 V / 50-60 Hz  |
| Power Consumption            | 50-80 VA  |
| Dimensions (WxLxH)           | 9.5 x 15.3 x 12.8 in / 241 x 390 x 324 mm   |
| Weight                       | 12.1 lb / 5.5 kg  |
| Clinical Test                | ISO81060-2: 2018  |
| Operating Conditions         | 50 to 104°F (10 to 40°C) / 15% to 85% RH / 70 to 106 kPa  |
| Storage/Transport Conditions | -4 to 140°F (-20 to 60°C) / 10% to 95% RH / 70 to 106 kPa   |

## Packaging and Shipping Specifications

|  | Dimensions (WxLxH)                      | Weight       | Case Cube         | UPC / GTIN   |
|--|---|--------------|-------------------|--|
| TM-2657WP<br>TM-2657WP-RS<br>TM-2657WP-BT<br>TM-2657WP-BLE | 22.8 x 18.1 x 13.3 in (59 x 48 x 36 cm) | 20 lb (9 kg) | 3.2 ft³ / 0.09 m³ | 0 93764 60430 0<br>0 93764 60464 5<br>0 93764 60465 2<br>0 93764 60466 9 |

## Optional Accessories

| Part Number                            | Order Code     |
|--|----------------|
| Arm Cuff Cover (Pack of 5)             | AX-134005759-S |
| Dedicated Stand                        | TM-ST520       |
| Adjustable Height Stool                | TM-STA001      |
| Printer Paper (Pack of 5)              | AX-PP147-S     |
| Disposable Kiosk Sleeves (Pack of 400) | AS-134010367   |

## What is in the Box

- 1 Kiosk monitor
- 1 Power cable
- 1 Printer paper roll
- 3 Instruction manuals
- 2 Patient instruction panels